Thank you for choosing this product.
Please read this manual carefully before using it. Please keep it in good condition for future inquiries.

Parameter description

Product name: Cervical vertebra massager
Rated power: 3.7w
Charging voltage: DC 5.0V

Common problems and handling

Serial number	phenomenon	possible reason	Approach
1	Massager can't start	Exhausted during standby	It can be started normally after charging.
2	The massager stops working	The timing message program has been completed or the battery power has been used up	Re-open the battery switch

If there is still an abnormality in the above processing, please contact your local dealer for repair service.

Use the power plug and power cord


1. Please check that the voltage does meet the product specifications.
2. Do not use other similar accessories instead of this product accessories.
3. Do not touch the power supply or plug with wet hands to avoid malfunction or electric shock.
4. Do not charge the product when the voltage is unstable.

Do not use this product when the following conditons occur


1. The massager is abnotmal, and the massage device is damaged or damaged.
2. If there is any serious or severe pain durion the massage, stop the massage immediately and seek medical advice.
3. In case of the above situation, if it is confirmed that the product is the cause, please send the product to the repair service center for inspection and maintenance.

Product configuration

Your U-shaped pillow includes: main unit, charging cable, instruction manual, certificate, warranty card

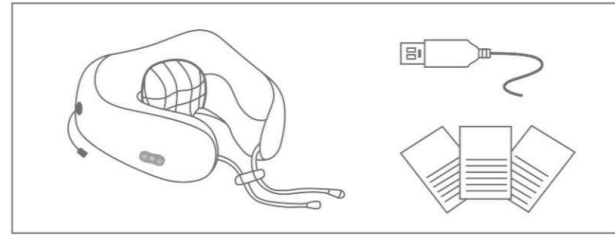Kneading type

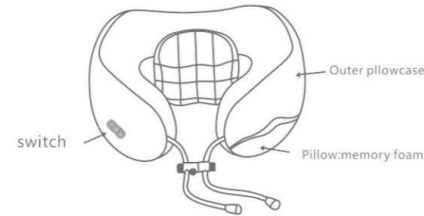Rope buckle design

Strong bond and weak bond

Positive Inversion

One button start

warranty service

Dear customer, thank you for your support. The products you purchased are produced under strict quality control and standard inspection process. This warranty card refers to the failure of the product under normal use conditions, according to the warranty card. Provisions to provide repairs and services to customers

Warranty description

1. From the date of purchase, if there is a malfunction under normal use, you can enjoy one year of free warranty service, except for damage (such as: improperly disassembling parts, handling damage, etc.), natural damage (such as: fire, earthquake, flood, Damage caused by rodent damage, etc. is not covered by this warranty.
2. If the warranty card is lost and will not be reissued, please keep it in a safe place.
3. No free warranty service due to malfunction caused by improper use of consumers themselves
4. If the product has expired within the warranty period, if the company still has accessories, it will be given warranty service, and certain parts and handling fees will be charged.
5. This warranty card is valid only in the territory of the Chinese Republic

The following conditions are not covered by the warranty

1. Failure to use the prohibited operation in accordance with the instructions
2. No purchase invoice or receipt, no valid warranty certificate, fake and sboddy products
3. Return the product packaging is incomplete, the accessories are not complete, or the user disassembles, repairs themselves, etc.
4. Losses caused by human factors, such as rubbing, immersion, violent impact, etc.
5. The product affects the seller's secondary sales due to human reasons
6. The information written on the warranty card or the warranty card does not match the actual product

U-shaped massage pillow

Instructions for the Use of Products

We reserve the right to improve our products and specifications without notice.

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